Anti-Vimentin Antibody [PDH0-01] IRS044



Hangzhou Huaan Biotechnology Co., Ltd.

Orders:0086-571-88062880

Technical:0086-571-89986345

5 Service mail:support@huabio.cn



11.

Applications:WB=Western blot IHC-P=Immunohistochemistry (paraffin) IF-Cell=Immunofluorescence (Cell) IF-Tissue=Immunofluorescence (Tissue) FC=Flow cytometry IP=Immunoprecipitation

IRS044 - Page 2

Images

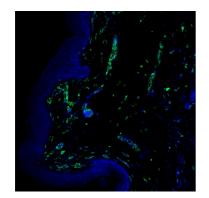


Fig1: mIHC analysis of human skin tissue (Formalin/PFA-fixed paraffin-embedded sections) with Rabbit anti-Vimentin antibody (IRS044) at 1/100 dilution. The immunostaining was performed with the IRISKit® HyperView mTSA Kit (MH900206). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0) for 30 mins at 95°C. DAPI (blue) was used as a nuclear counter stain. Image acquisition was performed with Olympus VS200 Slide Scanner.

Note: All products are "FOR RESEARCH USE ONLY AND ARE NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE".

Background References

- 1. Ridge KM et al. Roles of vimentin in health and disease. Genes Dev. 2022 Apr
- 2. Kuburich NA et al. Vimentin and cytokeratin: Good alone, bad together. Semin Cancer Biol. 2022 Nov

Hangzhou Huaan Biotechnology Co., Ltd.



Orders:0086-571-88062880

Technical:0086-571-89986345

Service mail:support@huabio.cn

Applications:WB=Western blot IHC-P=Immunohistochemistry (paraffin) IF-Celt=Immunofluorescence (Celt) IF-Tissue=Immunofluorescence (Tissue) FC=Flow cytometry IP=Immunoprecipitation