

# Human IL-8 ELISA Instructions

**Cat: EH0014**

## Content

	CAT	Volume
① CP (Coated Plate)	EH0014CP	96 well
② S (Standard)	EH0014S,S1-S7,S0	9 vial
③ DA-H (Detect Antibody-HRP)	EH0014DA-H	6 ml/bottle
④ AB (Assay Buffer 1×)	EAB01	12 ml/bottle
⑤ TS (TMB Substrate)	ETS01	12 ml/bottle
⑥ SS (Stop Solution)	ESS01	12 ml/bottle
⑦ WB (Wash Buffer 10×)	EWB01	50 ml/bottle
⑧ SF (Sealer Film)	ESF01	6 piece

**NOTE:** After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.

## REAGENT PREPARATION

### Washing Buffer (1×) Preparation


Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

### Standard Curve Preparation:

S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The Human IL-8 Standard EH0014S 30000 pg/ml 30  $\mu$ l + 270  $\mu$ l SPB serves as the high standard (3000 pg/ml). Pipette 200  $\mu$ l of SPB into each tube. Use the high standard to produce a 1:2 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

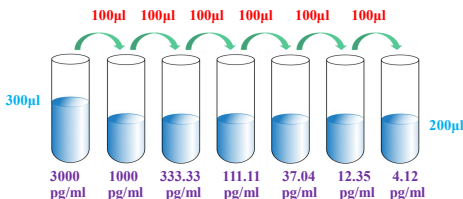


## ASSAY PROCEDURE

**Bring all reagents and samples to room temperature before use.**

- ① Prepare all reagents and working standards as directed in the previous sections.
- ② Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.
- ③ Add 50  $\mu$ l of **AB** (Assay Buffer) to each well.
- ④ Add 50  $\mu$ l or 10  $\mu$ l of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.
- ⑤ Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hours on a microplate **shaker** set at 500 rpm.
- ⑥ Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300  $\mu$ l). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.
- ⑦ Add 50  $\mu$ l of **DA-H** (Detect Antibody) to each well.
- ⑧ Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate **shaker** set at 500 rpm.
- ⑨ Repeat aspiration/**wash** as in step 6.
- ⑩ Add 100  $\mu$ l of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.
- ⑪ Add 100  $\mu$ l of **SS** (Stop Solution) to each well.
- ⑫ Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

## TYPICAL DATA

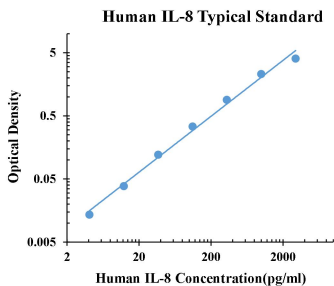

pg/ml	O.D.	Average	Corrected
0.00	0.0381	0.0373	0.0377
4.12	0.0493	0.0533	0.0513
12.35	0.0743	0.0785	0.0764
37.04	0.1583	0.1604	0.1594
111.11	0.3773	0.3767	0.3770
333.33	0.9294	0.9447	0.9371
1000.00	2.3256	2.3424	2.3340
3000.00	4.0171	4.1382	4.0777

## SENSITIVITY

The minimum detectable dose (MDD) of human IL-8 is typically less than 1.45 pg/ml (50  $\mu$ l of sample volume) or 3.63 pg/ml (10  $\mu$ l of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

## PRECISION

**Intra-assay Precision (Precision within an assay)** Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

**Inter-assay Precision (Precision between assays)**

Sample Number	Intra-assay Precision			Inter-assay Precision		
	S1	S2	S3	S1	S2	S3
22	22	22	22	6	6	6
Average (pg/ml)	54.3	279.3	896.1	51.2	259.7	780.4
Standard deviation	1.9	8.3	31.3	3.4	14.2	25.7
Coefficient of variation (%)	3.4	3.0	3.5	6.7	5.5	3.3

## RECOVERY

The spike recovery was evaluated by spiking 3 levels of human IL-8 into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 84% to 117% with an overall mean recovery of 102%.

## LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of IL-8 in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 94% to 111% with an overall mean recovery of 99%.

## SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy volunteers were evaluated for the presence of IL-8 in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	43.4-359.3	100	127.2

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.

## CALIBRATION

The NIBSC/WHO British Standard for human leukocyte IL-8 89/520 was evaluated in this kit. To convert sample values obtained with the Human IL-8 kit to relative approximate NIBSC units, use the equation below:

NIBSC/WHO (89/520) approximate value (U/ml) =  $0.001 \times$  Human IL-8 value (pg/ml)