

Human CD30 ELISA Instructions

Content

	CAT	Volume
 CP (Coated Plate) 	EH0105CP	96 well
Ø S (Standard)	EH0105S,S1~S7,S0	9 vial
3 DA (Detect Antibody)	EH0105DA	6 ml/bottle
SH (Streptavidin-HRP)	ESH01	12 ml/bottle
6 AB (Assay Buffer 1×)	EAB01	12 ml/bottle
6 TS (TMB Substrate)	ETS01	12 ml/bottle
SS (Stop Solution)	ESS01	12 ml/bottle
8 WB (Wash Buffer 10×)	EWB01	50 ml/bottle
9 SF (Sealer Film)	ESF01	6 pieces

NOTE: After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.

REAGENT PREPARATION

Washing Buffer (1×) Preparation


Pour entire contents (50 ml) of the Washing Buffer Concentrate (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

Standard Curve Preparation:

S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The human CD30 Standard EH0105S 40000 pg/ml 30 μ l + 270 μ l SPB serves as the high standard (4000 pg/ml). Pipette 200 μ l of SPB into each tube. Use the high standard to produce a 1:2 dilution series. Mix each tube thoroughly before the next transfer. SPB serves as the zero standard (0 pg/ml).

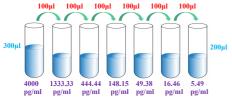

Cat:EH0105

ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

 Prepare all reagents and working standards as directed in the previous sections.

Remove excess CP (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.

3 Add 50 μl of AB (Assay Buffer) to each well.

4 Add 50 μl or 10 μl of Standard or sample per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.

S Add 50 μl of **DA** (Detect Antibody) to each well.

6 Cover with an SF (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hour on a microplate shaker set at 500 rpm.

7 Aspirate each well and wash, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 μ l). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and blot it against clean paper towels.

8 Add 100 μl of SH (Streptavidin-HRP) to each well.

Over with a new SF (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate shaker set at 500 rpm.

10 Repeat aspiration/wash as in step 7.


Add 100 μl of TS (TMB Substrate) to each well.
 Incubate for 5-30 minutes at room temperature.

12 Add 100 μl of SS (Stop Solution) to each well.

B Determine the optical density within 30 minutes, using microplate reader set to 450 nm corrected with 570 nm or 630 nm.



TYPICAL DATA

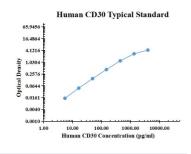

pg/ml	0.	O.D.		Corrected
0.00	0.0199	0.0202	0.0201	
5.49	0.0348	0.0353	0.0351	0.0150
16.46	0.0703	0.0685	0.0694	0.0494
49.38	0.1718	0.1691	0.1705	0.1504
148.15	0.4668	0.4436	0.4552	0.4352
444.44	1.2290	1.2290	1.2290	1.2090
1333.33	2.7640	2.8060	2.7850	2.7650
4000.00	4.3046	4.2674	4.2860	4.2660

SENSITIVITY

The minimum detectable dose (MDD) of human CD30 is typically less than 0.80 pg/ml (50 μ l of sample volume) or 3.19 pg/ml (10 μ l of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays)

	Intra-assay Precision		Inter-assay Precision		sion	
Sample	S1	82	S 3	S1	S2	S3
Number	22	22	22	6	6	6
Average (pg/ml)	87.5	398.2	1088.4	84.1	408.2	1263.3
Standard Deviation	5.2	21.2	81.8	2.9	25.6	82.2
Coefficient Of Variation (%)	6.0	5.3	6.0	3.5	6.3	6.5

RECOVERY

The spike recovery was evaluated by spiking 3 levels of human CD30 into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 83% to 98% with an overall mean recovery of 93%.

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of CD30 in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 89% to 97% with an overall mean recovery of 95%.

SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy volunteers were evaluated for the presence of human CD30 in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (pg/ml)	Detectable %	Mean of Detectable (pg/ml)
Serum	30	643.28-15 50.80	100	965.76

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.