

Human Complement C5 ELISA Instructions

Cat:EH0108

Content

	CAT	Volume
① CP (Coated Plate)	EH0108CP	96 well
② S (Standard)	EH0108S,S1-S7,S0	9 vial
③ SD (Sample Diluent)	ESD01	15ml/bottle
④ DA (Detect Antibody)	EH0108DA	6 ml/bottle
⑤ SH (Streptavidin-HRP)	ESH01	12 ml/bottle
⑥ AB (Assay Buffer 1×)	EAB01	12 ml/bottle
⑦ TS (TMB Substrate)	ETS01	12 ml/bottle
⑧ SS (Stop Solution)	ESS01	12 ml/bottle
⑨ WB (Wash Buffer 10×)	EWB01	50 ml/bottle
⑩ SF (Sealer Film)	ESF01	6 pieces

NOTE: After the kit is opened, the stabilization period of each content is 30 days, so please use the kit within 30 days after opening.

Sample Dilution

Samples such as serum 、 plasma require at least a 20000-fold dilution into Sample Diluent.

REAGENT PREPARATION

Washing Buffer (1×) Preparation


Pour entire contents (50 ml) of the **Washing Buffer Concentrate** (10×) into a clean 500 ml graduated cylinder. Bring to final volume of 500 ml with glass-distilled or deionized water. Transfer to a clean wash bottle and store at 2 to 25°C.

Standard Curve Preparation:

S1 to S7 and S0 is ready to use for serum and plasma.

Other sample type, prepare the standard curve with whatever buffer (SPB, Sample Prepared Buffer) is used to prepare the sample, such as cell culture supernatant, tissue grinding liquid, cell lysate, etc. Urine sample use AB (Assay Buffer) prepare standard curve.

The human Complement C5 Standard EH0108S 20000 pg/ml 30 μ l + 270 μ l SPB serves as the high standard (2000 pg/ml). Pipette 200 μ l of Sample Diluent into each tube. Use the high standard to produce a 1:2 dilution series. Mix each tube thoroughly before the next transfer. Sample Diluent serves as the zero standard (0 pg/ml).

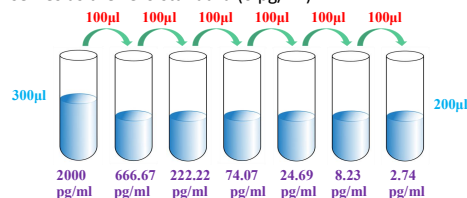

ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

① Prepare all reagents and working standards as directed in the previous sections.

② Remove excess **CP** (Coated Plate) strips from the plate frame, return them to the foil pouch and reseal.

③ Add 50 μ l of **AB** (Assay Buffer) to each well.

④ Add 50 μ l or 10 μ l of **Standard or sample** per well. Ensure reagent addition is uninterrupted and completed within 15 minutes.

⑤ Add 50 μ l of **DA** (Detect Antibody) to each well.

⑥ Cover with an **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 1 hour on a microplate **shaker** set at 500 rpm.

⑦ Aspirate each well and **wash**, repeating the process four times. Wash by filling each well with WB (Washing Buffer 300 μ l). Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining WB (Washing Buffer) by aspirating or decanting. Invert the plate and **blot** it against clean paper towels.

⑧ Add 100 μ l of **SH** (Streptavidin-HRP) to each well.

⑨ Cover with a new **SF** (Sealer Film). Incubate at room temperature (18 to 25°C) for 30 min on a microplate **shaker** set at 500 rpm.


⑩ Repeat aspiration/**wash** as in step 7.

⑪ Add 100 μ l of **TS** (TMB Substrate) to each well. Incubate for 5-30 minutes at room temperature.

⑫ Add 100 μ l of **SS** (Stop Solution) to each well.

⑬ Determine the optical density within 30 minutes, using microplate **reader** set to 450 nm corrected with 570 nm or 630 nm.

TYPICAL DATA

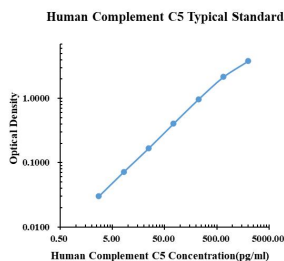

pg/ml	O.D.	Average	Corrected
0.00	0.0191	0.0181	0.0186
2.74	0.0494	0.0480	0.0487
8.23	0.0899	0.0917	0.0908
24.69	0.1898	0.1834	0.1866
74.07	0.4362	0.4138	0.4250
222.22	0.9751	0.9995	0.9873
666.67	2.2140	2.1390	2.1765
2000.00	3.8600	3.8120	3.8360

SENSITIVITY

The minimum detectable dose (MDD) of human Complement C5 is typically less than 0.05 pg/ml (50 µl of sample volume) or 3.14 pg/ml (10 µl of sample volume).

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

Intra-assay Precision (Precision within an assay) Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays)

	Intra-assay Precision			Inter-assay Precision		
Sample Number	S1	S2	S3	S1	S2	S3
	22	22	22	6	6	6
Average (pg/ml)	17.8	92.6	287.7	16.6	85.8	261.2
Standard Deviation	0.7	4.5	15.2	0.7	4.8	16.9
Coefficient of Variation (%)	3.7	4.9	5.3	4.5	5.6	6.5

RECOVERY

The spike recovery was evaluated by spiking 3 levels of human Complement C5 into health human serum sample. The un-spiked serum was used as blank in this experiment.

The recovery ranged from 80% to 101% with an overall mean recovery of 92%.

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of Complement C5 in human serum and diluted with Sample Diluent to produce samples with values within the dynamic range of the assay.

The linearity ranged from 83% to 95% with an overall mean recovery of 90%.

SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy volunteers were evaluated for the presence of human Complement C5 in this assay. No medical histories were available for the donors.

Sample Matrix	Sample Evaluated	Range (µg/ml)	Detectable %	Mean of Detectable (µg/ml)
Serum	30	1.79-58.93	100	25.41

n.d. = non-detectable. Samples measured below the sensitivity are considered to be non-detectable.