Anti-Granzyme B Antibody

IRS054



Product Type: Rabbit polyclonal IgG, primary antibodies

Species reactivity: Mouse
Applications: mIHC

Molecular Wt: Predicted band size: 28 kDa

Description: Abundant protease in the cytosolic granules of cytotoxic T-cells and NK-cells which activates

caspase-independent pyroptosis when delivered into the target cell through the immunological synapse. It cleaves after Asp. Once delivered into the target cell, acts by catalyzing cleavage of gasdermin-E (GSDME), releasing the pore-forming moiety of GSDME, thereby triggering pyroptosis and target cell death. Seems to be linked to an activation cascade of caspases (aspartate-specific cysteine proteases) responsible for apoptosis execution. Cleaves caspase-3, -7, -9 and 10 to give rise to active enzymes mediating

apoptosis.

Immunogen: Recombinant protein within human Granzyme B aa 50-247.

Positive control: Mouse osteosarcoma tissue, mouse spleen tissue.

Subcellular location: Secreted.

Database links: SwissProt: P04187 Mouse

Recommended Dilutions:

mI HC 1:100


Storage Buffer: PBS (pH7.4), 0.1% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide.

Storage Instruction: Store at +4℃ after thawing. Aliquot store at -20℃ or -80℃. Avoid repeated freeze / thaw

cycles.

Purity: Immunogen affinity purified.

Hangzhou Huaan Biotechnology Co., Ltd.

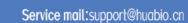



Images

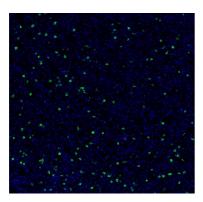


Fig1: mIHC analysis of mouse osteosarcoma tissue (Formalin/PFA-fixed paraffin-embedded sections) with Rabbit anti-Granzyme B antibody (IRS054) at 1/100 dilution. The immunostaining was performed with the IRISKit® HyperView mTSA Kit (MH900206). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0) for 30 mins at 95 $^{\circ}$ C. DAPI (blue) was used as a nuclear counter stain. Image acquisition was performed with Olympus VS200 Slide Scanner.

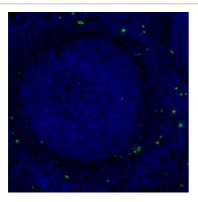


Fig2: mIHC analysis of mouse spleen tissue (Formalin/PFA-fixed paraffin-embedded sections) with Rabbit anti-Granzyme B antibody (IRS054) at 1/100 dilution. The immunostaining was performed with the IRISKit® HyperView mTSA Kit (MH900206). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0) for 30 mins at $95\,^{\circ}$ C. DAPI (blue) was used as a nuclear counter stain. Image acquisition was performed with Olympus VS200 Slide Scanner.

Note: All products are "FOR RESEARCH USE ONLY AND ARE NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE".

Background References

1. Hameed A. et. al. Characterization of three serine esterases isolated from human IL-2 activated killer cells. J. Immunol. 141:3142-3147(1988).