Anti-MAP2 Antibody

IRS071



Product Type: Rabbit polyclonal IgG, primary antibodies

Species reactivity: Mouse
Applications: mIHC

Molecular Wt: Predicted band size: 200 kDa

Description: This gene encodes a protein that belongs to the microtubule-associated protein family. The

proteins of this family are thought to be involved in microtubule assembly, which is an essential step in neurogenesis. The products of similar genes in rat and mouse are neuron-specific cytoskeletal proteins that are enriched in dentrites, implicating a role in determining and stabilizing dentritic shape during neuron development. A number of alternatively spliced

variants encoding distinct isoforms have been described.

Immunogen: Synthetic peptide within human MAP2 aa 1760-1810.

Positive control: Mouse brain tissue.

Subcellular location: Cytoskeleton, dendrite.

Database links: SwissProt: P20357 Mouse

Recommended Dilutions:

mIHC 1:100

Storage Buffer: PBS (pH7.4), 0.1% BSA, 40% Glycerol. Preservative: 0.05% Sodium Azide.

Storage Instruction: Store at $+4^{\circ}$ C after thawing. Aliquot store at -20° C or -80° C. Avoid repeated freeze / thaw

cycles.

Purity: Immunogen affinity purified.

Hangzhou Huaan Biotechnology Co., Ltd.





Images

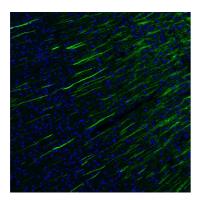


Fig1: mIHC analysis of mouse brain tissue (Formalin/PFA-fixed paraffin-embedded sections) with Rabbit anti-MAP2 antibody (IRS071) at 1/100 dilution. The immunostaining was performed with the IRISKit® HyperView mTSA Kit (MH900206). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0) for 30 mins at $95\,^{\circ}\mathrm{C}$. DAPI (blue) was used as a nuclear counter stain. Image acquisition was performed with Olympus VS200 Slide Scanner.

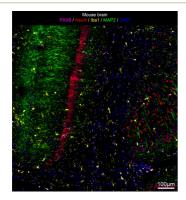


Fig2: mIHC analysis of mouse brain tissue (Formalin/PFA-fixed paraffin-embedded sections) with PAX6 (IRS066), NeuN, Iba1 (IRS068) and MAP2 (IRS071) antibody at 1/100 dilution. The immunostaining was performed with the IRISKit® HyperView mTSA Kit (MH900206). Heat mediated antigen retrieval with Tris-EDTA buffer (pH 9.0) for 30 mins at 95 °C. DAPI (blue) was used as a nuclear counter stain. Image acquisition was performed with Olympus VS200 Slide Scanner.

Note: All products are "FOR RESEARCH USE ONLY AND ARE NOT INTENDED FOR DIAGNOSTIC OR THERAPEUTIC USE".

Background References

- 1. DeGiosio R. et. al. MAP2 immunoreactivity deficit is conserved across the cerebral cortex within individuals with schizophrenia. NPJ Schizophr. 2019 Aug
- 2. Gumy LF. et. al. MAP2 Defines a Pre-axonal Filtering Zone to Regulate KIF1- versus KIF5-Dependent Cargo Transport in Sensory Neurons. Neuron. 2017 Apr